



HARVARD®

Harvard LiquiCord®

Retraction paste for temporary gingival displacement

Harvard LiquiCord® is a retraction paste with hemostatic effect for temporary gingival displacement. Depending on the clinical situation and the respective mode of operation, it can be used directly or with conventional retraction cord. If necessary with the aid of compression caps. Harvard LiquiCord® is well compatible with vinylpolysiloxane impression materials.


For the Dentist:

- Easy to use
- Safe application
- Time saving

For the Patient:

- Gentle to the mucosa
- Pleasant
- Neutral taste

A supplement or an alternative to thread ...

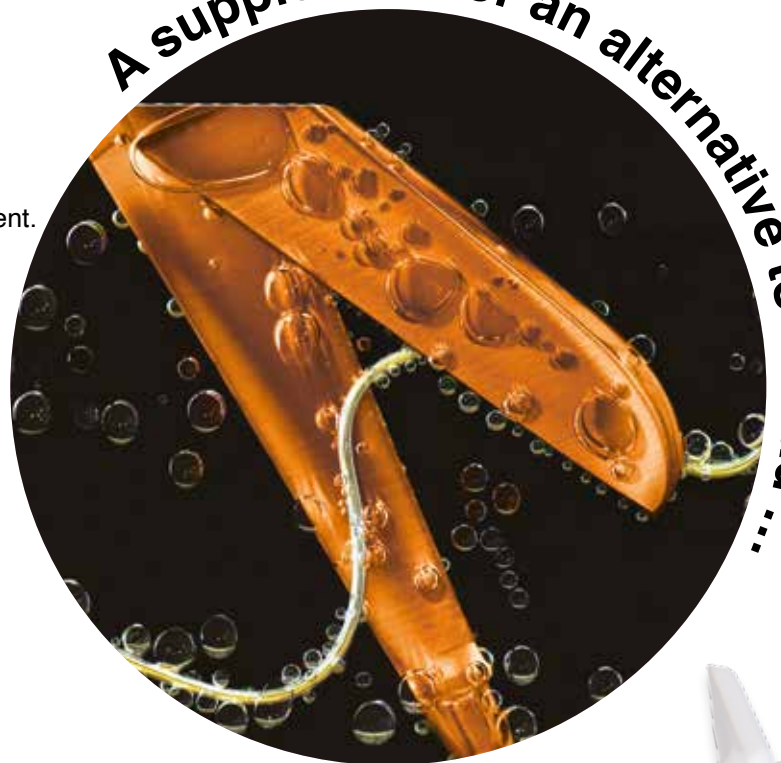

Optional accessories: Suitable Sulcus Former



small



medium



large




Application of Harvard LiquiCord® ...

... for accurate impressions the gingiva retraction is very important: after a tooth preparation the gingival sulcus normally is filled with saliva, blood or other liquids and therefore not accessible for the impression material. This requires effective retraction with Harvard LiquiCord®. It temporarily displaces the gingiva and keeps the sulcus dry and protected. In addition Harvard LiquiCord® has a hemostatic effect.

Harvard LiquiCord® can be used for all indications of temporary displacement of the marginal gingiva and to provide a dry and expanded sulcus:

- For digital impressions
- For conventional impressions
- Impressions for restorations (temporary and permanent)
- Preparations of class II and V

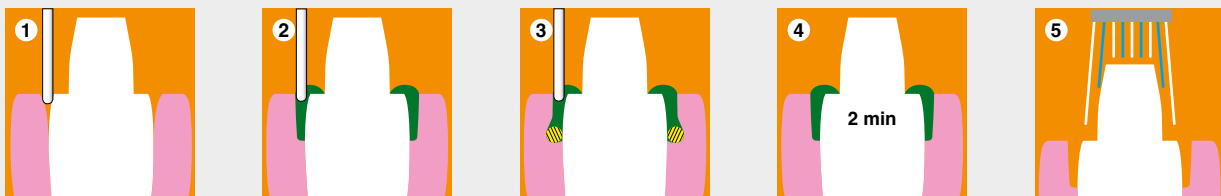

1. Introduce the tip/needle into the sulcus, opening the sulcus.
2. Move the top of the tip slowly and evenly around the tooth while pressing out the paste. Fill the sulcus all around with enough Harvard LiquiCord so that an excess appears.
3. If necessary, Harvard LiquiCord can be used in combination with a retraction cord, that is placed before in the sulcus, to open the sulcus further.
4. Let Harvard LiquiCord take effect for at least 2 minutes and keep away all moisture from the site during this time.
5. After the residence time remove the Harvard LiquiCord completely from the sulcus with air and water and with the aid of a suction device.



Order information	Order no.
Harvard LiquiCord – syringe 10 x 0.7 g, 20 needle tips	7083870
Harvard Sulcus Former – small Refill bag with 100 Sulcus Former	7083872
Harvard Sulcus Former – medium Refill bag with 100 Sulcus Former	7083873
Harvard Sulcus Former – large Refill bag with 50 Sulcus Former	7083874
Harvard NeedleTips H18 white Refill bag with 50 needle tips	7083875

1892
seit
HARVARD



HARVARD®

Marke und Qualität seit 1892

Harvard Dental International GmbH
Margaretenstr. 2 - 4, 15366 Hoppegarten, Germany
Phone: + 49 (0) 30/99 28 978-0
Fax: + 49 (0) 30/99 28 978-19
info@harvard-dental-international.de

Liability is excluded for all printing errors and omissions. Before using our Harvard products, the respective directions for use should be noticed in every case. All measurements are internal measurements of Harvard Dental International.

Harvard Distribution Partner.